



A NEW WAY  
TO HELP YOUR  
PATIENT FEEL SAFE

Navina Fecal Incontinence Insert is an innovation  
designed to prevent leakage accidents

**Navina**™


# SHIFTING FROM WORRY TO NEW HOPE

What is quality of life, really? Ask someone with bowel issues who is prevented from leaving the house due to fear of leakage and embarrassing accidents. When everything revolves around staying close to a toilet, life becomes very limited. In fact, it can lead to delayed treatment, self-imposed isolation and depression.

Wellspect presents a smart solution designed to help manage the symptoms of fecal incontinence. This new innovation is the soft and flexible Navina Fecal Incontinence Insert, made of silicone and liquid. It uses the science of fluid dynamics to give the insert an anatomic fit and help prevent involuntary loss of stool.

“I don’t say this very often, but I’m so, so very happy with this product. I can really see it making a real difference in my patients’ lives.”


Bev Collins, Clinical Nurse Lead,  
Wellspect HealthCare, UK

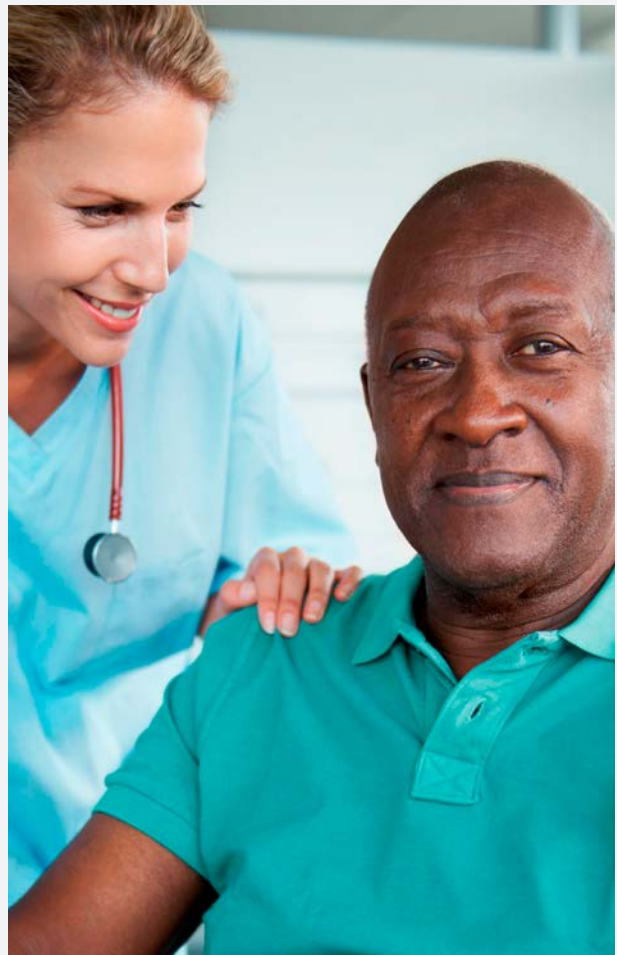

## STEP BY STEP INSTRUCTION



**1.** Hold the pre-lubricated insert by the white plastic applicator and align the tip of the bulb with the anus.

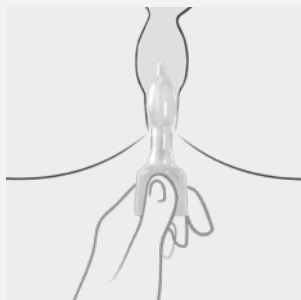

**2.** Gently push the insert until the broader retainer rests on the anus. Withdraw the applicator and discard it in the regular trash.



**3.** The liquid-filled insert adapts to the contours of the body and helps prevent leakage.



**4.** To remove it, grasp the retainer and gently pull to remove the insert and discard it with the regular trash.



## AN INNOVATIVE DESIGN FOR MORE CONFIDENCE

The flexible, liquid-filled insert is engineered to conform, move and adapt as the fluid is transferred throughout the device.

This makes it easy to insert and remove, and comfortable to wear. Navina Fecal Incontinence Insert is designed to make it easy to handle, prevent leakage and allow for an independent and active life.

- Engineered to Protect against accidental bowel leakage
- Designed to be Comfortable
- Intended to be Hygienic

### 24 HOUR USE

It may be worn continuously for up to 24 hours

### HYGIENE AND EASE-OF-USE

The plastic applicator keeps the insert steady during insertion and helps keep fingers clean.

### FLEXIBILITY AND COMFORT

Made from soft silicone and filled with smooth liquid, the insert is carefully crafted to ensure comfortable insertion and to adapt and conform to movement.

### SECURE POSITIONING

The broad external retainer collapses comfortably and is designed to hold the device securely in place.

### DUAL SEAL

First, the bulb rests against the anorectal junction, providing a seal to prevent accidental fecal discharge. Second, the shaft conforms to and fills the anal canal allowing pressure from the sphincters to seal the canal in order to prevent bypass leakage.

**GET IN TOUCH TODAY.  
YOUR PATIENTS WILL  
BENEFIT TOMORROW.**

If you have questions or need patient brochures, please contact Wellspect Customer Service.  
[www.wellspect.com](http://www.wellspect.com).

## NAVINA FECAL INCONTINENCE INSERT - DESCRIPTION:

The Fecal Incontinence Insert is a soft, flexible, liquid filled insert designed to prevent fecal incontinence. The device is non-sterile and for single use only.

### INTENDED USE

The rectal insert is designed for self-insertion to seal and help prevent the involuntary leakage of stool from the rectum.

### INDICATIONS

The Fecal Incontinence Insert is indicated for the management of Accidental Bowel Leakage (ABL) due to bowel incontinence.

### CONTRAINDICATIONS

Do not use Navina Fecal Incontinence Insert in the pediatric population.

At Wellspect we develop innovative continence care solutions that improve quality of life for people with bladder and bowel problems. We inspire our users to build self-confidence and independence as well as good health and well-being. We have been leading the industry for over 40 years with our product brands LoFric<sup>®</sup> and Navina<sup>™</sup>. We always aim to minimize the environmental impact of our products and passionately strive to become climate neutral.

We work together with users and healthcare professionals to improve clinical outcome in a sustainable way, now and for the future.

**Wellspect. A Real Difference.**

For more information about our products and services, please visit [Wellspect.com](https://www.wellspect.com).

Join the conversation on LinkedIn, Twitter, Facebook and Instagram.

[wellspect.com](https://www.wellspect.com)



Minnesota Medical Technologies  
2446 Henry Road Northwest  
Stewartville, MN 55976 USA  
507-533-0366  
[www.mnmedtech.us](https://www.mnmedtech.us)



Importer: Wellspect HealthCare  
(A division of Dentsply IH AB)  
Aminogatan 1  
PO Box 14  
SE-431 21 Mölndal  
Sweden  
[www.wellspect.com](https://www.wellspect.com)



MPS Medical Product Service GmbH  
Borngasse 20  
Braunfels 35619  
Germany